


# Hardy **CHROM**<sup>TM</sup> **CRE**

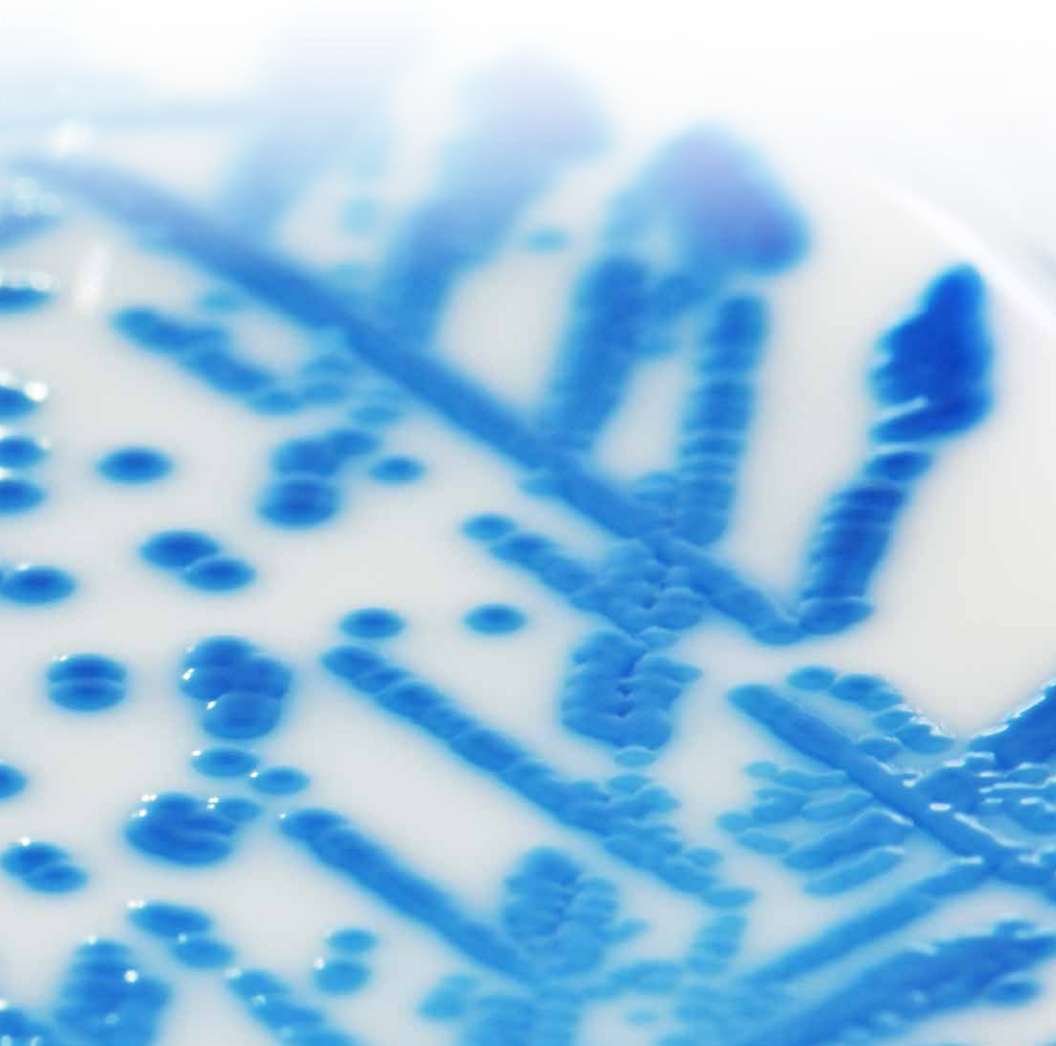


**Hardy Diagnostics is proud to introduce HardyCHROM CRE for the screening of *E. coli*, *K. pneumoniae*, *K. aerogenes*, *K. oxytoca*, *Enterobacter cloacae* complex, and *Serratia marcescens*.**

HardyCHROM™ CRE is a selective and differential chromogenic agar medium intended for the qualitative and presumptive detection from stool specimens of *Escherichia coli* that are non-susceptible to carbapenems as pink colonies and KES (*Klebsiella aerogenes*, *Klebsiella oxytoca*, *Klebsiella pneumoniae*, *Enterobacter cloacae* complex, and *Serratia marcescens*) that are non-susceptible to carbapenems as blue colonies.


HardyCHROM™ CRE is intended as an aid in the detection, identification of colonization and control of these bacteria in a healthcare setting. HardyCHROM™ CRE is not intended to diagnose infection or guide therapy. Results can be interpreted after incubation for 18-24 hours. Subculture to non-selective medium is required for confirming identification, antimicrobial susceptibility testing and epidemiological typing.

15x100mm, 10/pk, Cat. no.....G323



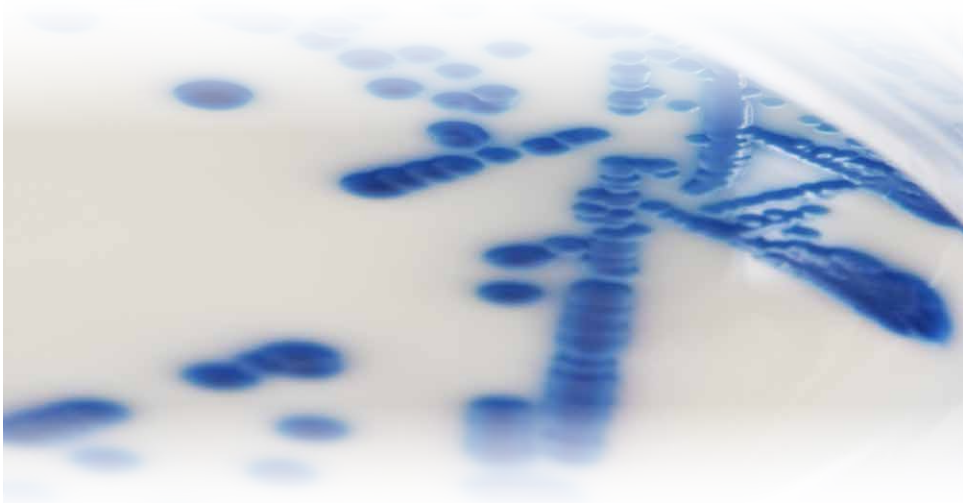

*Klebsiella*, *Enterobacter*, and *Serratia* spp. produce large, dark blue colonies (with or without pink halos) are presumptive positive for carbapenem non-susceptible KES (*Klebsiella aerogenes*, *Klebsiella oxytoca*, *Klebsiella pneumoniae*, *Enterobacter cloacae* complex, and *Serratia marcescens*).

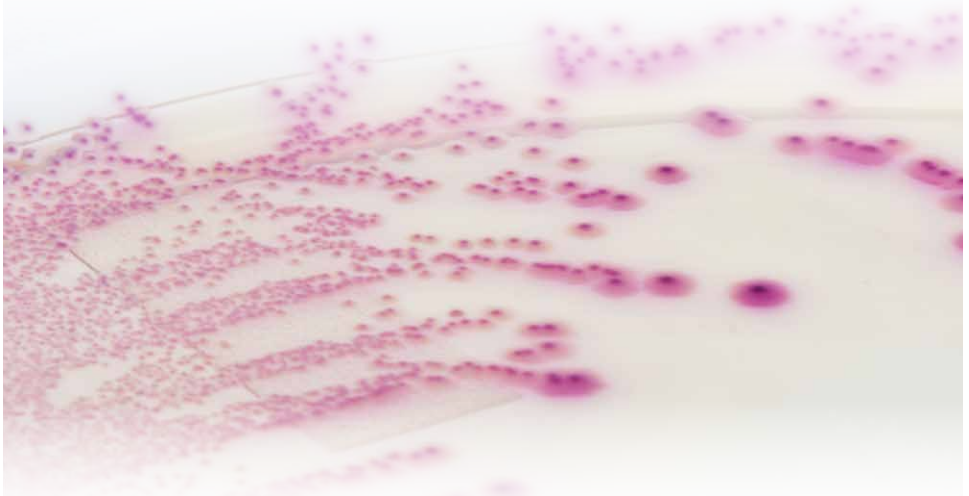

Non-susceptible *E. coli* produces colonies that are rose to magenta in color with darker pink centers.

Colonies that are not pink to magenta, blue, or blue with pink halos are negative. No carbapenem non-susceptible *Escherchia coli* or KES detected.

**Call 800.266.2222 or visit  
HardyDiagnostics.com for a  
free sample!**



**100%**  
Employee Owned



Made in the USA



FM 572526

Hardy Diagnostics has a Quality Management System that is certified to ISO 13485 and is a FDA licensed medical device manufacturer.

#### Headquarters

1430 West McCoy Lane  
Santa Maria, CA 93455  
800.266.2222  
Sales@HardyDiagnostics.com  
HardyDiagnostics.com

#### Distribution Centers

Santa Maria, California  
Olympia, Washington  
Salt Lake City, Utah  
Phoenix, Arizona  
Dallas, Texas  
Springboro, Ohio  
Lake City, Florida  
Albany, New York  
Raleigh, North Carolina

Copyright © 2019 Hardy Diagnostics

050319sg